

PRODUCT DATASHEET

Disposable Urine Bottle, Round, 1 L, w/Lid & Stick, WardPak™

Brand	: CURAS®
Category	: Nursing
Range	: Hygiene Pulp™
Product Group	: Disposable Urine Bottles
REF	: 08057

PRODUCT DESCRIPTION

Disposable round urine bottle with a water-resistant lid made of pulp material. 1 L capacity. Measuring sticks are included for emergency use. Suitable for disposal in a macerator. Packed in WardPak™, compact dispensing carton for easy dispensing, non-sterile.

INTENDED USE


Disposable round urine bottle for use by immobile male patients for the collection of urine. The bottle is also ideal for collecting urine drained from a urine bag. Suitable for single-use only.

FEATURES & BENEFITS

- Disposable and single-use can reduce the risk of cross-contamination and hospital-acquired infections.
- Water resistant, for a minimum period of 4 hours (applicable to urine bottle and lid only).
- Ergonomic jug design for convenient use and holding optimal volume.
- With a contoured handle for grip to ease handling during urination.
- Reinforced flat base provides extra strength and allows stable placement.
- Designed to easily fit into bedside urine bottle holder.
- Disposable lid helps to contain odours and prevent accidental spills.
- Disposable graduated stick for urine measurement purpose.
- All in one package for convenient one-stop access to the necessary related products.
- WardPak™, dispensing carton protects product from dust and contamination.
- Perforated opening on dispensing carton offers an easy dispensing option.
- Can be used with superabsorbent polymer (SAP) by Curas™ (REF 08100) to prevent accidental spills, control odour and improve waste handling.

PRODUCT SPECIFICATIONS

REF	Capacity	Dimensions (LxWxH)	Average Weight/Pc	Packing
08057	1 L	24.5x11x12 cm	Urine bottle: 50 ± 5 g Lid: 3-4 g Stick: ~2.5 g	WardPak™

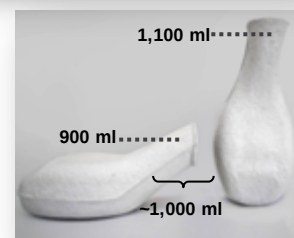

Composition

- Urine bottle and lid made of 100% recycled newsprint paper with water-based fluid resistant compound.
- Graduated stick made of cardboard.

Country of Origin

- China

Sterilization

- Non-sterile
- Can be sterilized by EtO and Gamma radiation.

Shelf Life

- 8 years

Selected Standards & Guidelines

- MDD 93/42/EEC (CE Marking, Class I), as amended by Directive 2007/47/EC.
- ISO 9001:2008
- ISO 2859-1 Sampling Plan
- PAS29:1999 Fluid resistance test for a minimum period of 4 hours (water at 35 ± 3 °C).
- PAS29:1999 Maceration test with a hospital sluice room type macerator, normal cycle.

Recommended Storage Condition

- Store at room temperature and avoid exposure to extreme temperatures.
- Keep at a cool dry place, away from moisture, heat and sources of ignition.

Proposed Disposal Method

- If regular waste, dispose of used product together with content into a pulp macerator, or in accordance with the local waste protocol. If a macerator is not available, incinerate or dispose the product into a landfill.
- If content is contaminated or contagious, treat as clinical waste and dispose it in accordance with the clinical waste protocol.

Product Symbols

PACKAGING INFO

REF	EAN-13 Barcode	WardPak™ (pcs/pack)	Carton (pcs/carton)	Gross Weight
08057	5051400080579	Urine Bottle: 28 Lids: 5 Graduated Stick: 10	Urine Bottle: 84 Lids: 15 Graduated Stick: 30	~2.1 kg

Components/Material Properties			Recommendations on Waste Handling				
Packaging	Material	Contains PVC	Incinerate	Landfill	Macerate	Recycle	Reuse
WardPak™	White Corrugated Carton	--	(√)	(√)		√	
Carton	Corrugated Carton					(√)	√
Transport Pallet	Wood	--				(√)	√

√ = Highly recommended

(√) = Recommended

ADDITIONAL INFO

- Curas® Hygiene Pulp™ products have been tested in several macerator models and are compliant under recommended operating procedures.
- HS Code: 4823709000
- Documentation available upon request. Kindly contact: info@curas.com