

VENOSAN® Eureka

Donning and doffing aid

January 2024


1. Composition of the product

VENOSAN® Eureka On consist of polyamide, cotton, PU coating and a PU pin

VENOSAN® Eureka Off consist of polyamide and a PU coating

REF	Product designation	Version
5002, 7312	VENOSAN® Eureka On	Normal/Plus
7143, 7313	VENOSAN® Eureka Off	Normal/ Plus

Sizing table

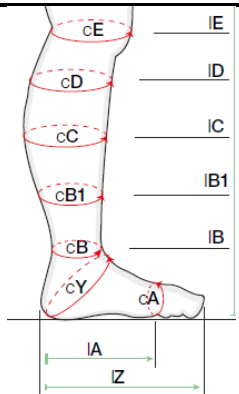VENOSAN® Eureka and Eureka Plus				
	cA	cY	length	
Eureka On	< 25cm	< 40cm	39,5cm	
Eureka On Plus	25-30cm	40-45cm	42cm	
Eureka Off	< 25cm	< 40cm	57cm	
Eureka Off Plus	25-30cm	40-45cm	78cm	

2. Packaging, structure and composition

2.1 Unit container

See 2.2

2.2 Shelf Container

Eureka On:

1 piece = a shelf container contains:

1 PE slider grip bag, IFU, Inlay and a PU pin.

Eureka Off:

1 piece= a shelf container contains:

1 PE slider grip bag, inlay and IFU.

According to the suppliers information, we currently have no information on ingredients of concern.

2.3 Transit Container

n/a

3. **Manufacture**

Eureka is manufactured in CZ.

4. **Description**

Eureka is intended to aid donning and doffing of all types of medical compression stockings.

5. **Properties**

Eureka is available in a donning (Eureka On) and Doffing (Eureka Off) aid. Eureka On and Off are in standard or Plus version.

6. **Intended purpose**

Eureka on and off can be used for easing donning and doffing MCH with both open and closed toe. Eureka mechanically expands the elastic material of the stocking in the area of the foot and lower leg to make it easier to pull on and take off the medical compression stockings.

7. **Medical device classification**

Eureka are medical devices of class I, Rule I according to MDR (EU) 2017/745.

8. **Biological evaluation and biocompatibility**

There is no known risk involved in the use of the medical device Eureka and it is designed, manufactured and packaged in such a way that it will not compromise the clinical condition or the safety of patients, or the safety and health of users and other persons when used under the conditions and for the purposes intended.

9. **Stability**

Store dry, dust free, at room temperature and protect from direct sunlight.
The material does not have a shelf life.


10. **Disposal**

Eureka can be disposed of in household waste. If local national regulations deviate, these instructions have to be followed.

11. **Legal manufacturer**

Swisslastic AG hereby confirms the accuracy of the information provided in this product data sheet.

12.01.2024

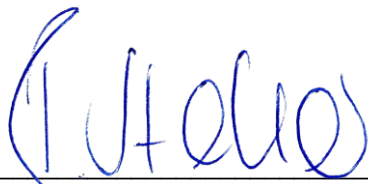

Date, Signature

Ramona Stecher

Product Management Swisslastic AG